
# HEIDELBERGER 30 / 75 / 140 / 190

<b>EN</b>	Instructions for use	<b>FR / B / CH</b>	Mode d'emploi	<b>RO</b>	Instrucțiuni de utilizare
<b>DE / AT / CH</b>	Gebrauchsanweisung	<b>HR</b>	Upute za uporabu	<b>SK</b>	Návod na použitie
<b>BG</b>	Инструкция за употреба	<b>HU</b>	Használati utasítás	<b>SL</b>	Navodila za uporabo
<b>CS</b>	Návod k použití	<b>IT</b>	Istruzioni per l'uso	<b>SV</b>	Bruksanvisning
<b>DK</b>	Bruksanvisning	<b>LT</b>	Naudojimo instrukcija	<b>NO</b>	Bruksanvisning
<b>GR / CY</b>	Οδηγίες χρήσης	<b>LV</b>	Lietošanas instrukcija		
<b>ES</b>	Instrucciones de uso	<b>NL / B</b>	Gebruikersinformatie		
<b>EE</b>	Kasutusjuhised	<b>PL</b>	Instrukcja użytkowania		
<b>FI</b>	Käyttöohje	<b>PT</b>	Instruções de utilização		

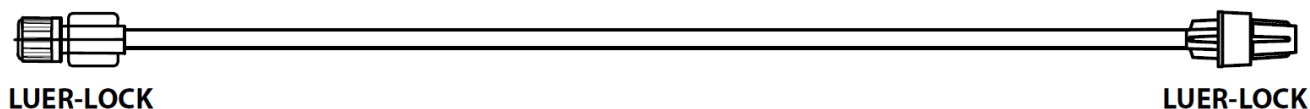REF:		
HEIDELBERGER 30	<b>CPHF51030DF</b>	EFB004AA030
HEIDELBERGER 75	<b>CPHF51070DF</b>	EFB004AA075
HEIDELBERGER 140	<b>CPHF51150DF</b>	EFB004AA140
HEIDELBERGER 190	<b>CPHF51300DF</b>	EFB004AA190

		HEIDELBERGER [mm]
ID / Вътрешен диаметър / Vnitřní průměr / Εσωτερική διάμετρος / Diámetro interior / Sisäläbimõõt Sisähalkaisija / DI / unutarnji promjer Belső átmérő / Vid. Skersm / Iekšējais diametrs / Binnendiam. Identyfikator / Vnútný priemer / NP		3
OD / AD / Външен диаметър / Vnější průměr / UD / Εξωτερική διάμετρος / Diámetro exterior / Välisläbimõõt Ulkohalkaisija / DE / vanjski promjer / Külső átmérő / Išor. skersm. / Ārējais diametrs / Buitendiam. / Średnica zewnętrzna / Vonkajší priemer / ZP		4 (4.1)

## EN / IRL / MT

**Important notice:** Please read the information on the primary packaging. Do not use if package is damaged. Sterilized by EO.

**Designated purpose:** The Heidelberg extension tube is designed to connect an infusion set and an intravenous cannula, or other infusion lines for Infusion Therapy.

**Indications:** Heidelberg extension tube for connection to infusion and transfusion sets. The Heidelberg extension tube is connected between the infusion set tubing and the intravenous cannula.

### Instructions for use:

1. while observing hygienic conditions, remove the product from the packaging. 2. Remove the cover of the positive LL connector.
3. Connect the LLM to the LLF connector of the set.
4. Connect the LLM connector to the intravenous cannula, or another connector of the infusion line.

**Patient population:** The product can be used in all patients who are prescribed infusion therapy. It has no gender or age restrictions.

**Designated user:** The device is intended for use by authorized persons, eg: health care professionals, licensed and/or board-certified physicians, medical professionals, and/or specially trained patients or caregivers (according to local regulations).

**Contraindications:** Medicines / solutions for which incompatibility with any of the materials used has been proven. The Summary of Product Characteristics of the respective drugs/solutions should be taken into account.

**Residual risks / Side effects:** General risks of infusion therapy including air embolism, local or systemic infection, contamination by foreign particles. Another hazard to patients, visitors or staff may be the leakage of cytotoxic or biological fluids.

**Safety measures:** Visual inspection of Heidelberg tubing before use for damage. – Do not use if protective caps are missing or loose. When using a positive pressure infusion, first remove the air from the container. Follow general safety precautions. Follow directions for use for parallel infusions. Please also take into account the information in the summary of product characteristics of the respective drugs/solutions regarding possible incompatibilities of these drugs/solutions. Check the tightness of all connections throughout the infusion.

**Warning:** Reuse of single-use devices poses a potential risk to the patient or user. It can lead to contamination or deterioration of the device's functionality, which can cause injury, illness, or death of the patient. Administration of light-sensitive drugs requires the use of UV protected intravenous administration sets.

**Period of use:** The duration of use depends on the treatment of the patient according to the Summary of Product Characteristics of the drug/solution. Replace the application set according to national guidelines.

**Disposal:** Disposal according to local guidelines and/or clinical protocols.

**Note for users:** Any serious incident that occurs in connection with this device must be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient resides.

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